

MOVIX 4/8 DREAM

MOVIX 4/8 DREAM



SEE THOROUGHLY, ENSURE SAFETY

MOVIX 4/8 DREAM

SYSTEM OVERVIEW

MOVIX 4/8 DREAM is a lightweight mobile x-ray unit (only 94 kg) equipped with a high frequency 4 kW or 8 kW generator (125 kV, 100 mA, 250 mAs) and with single or dual Portable Flat Panel Detector (FPD) that provides high image quality in any situation. An integrated computer enables **instantaneous** image display, immediate image processing, local and USB device storage, networked archiving and printing.

The lightweight portable FPD is as easy to use as a conventional X-ray cassette. When procedures are completed, the portable FPD is simply placed back the convenient storage compartment. **MOVIX 4/8 DREAM** is designed either for outdoor procedures or for in-room examinations.



COMPACTNESS

The dimensions of the **MOVIX 4/8 DREAM** (715 x 670 x 1410 mm) facilitate the transport in a vehicle or packed into its crate. On site, the mobile can be deployed in 3 steps, ready for exposures in less than 2 minutes, and provides image preview with 3 seconds.

Its manoeuvrability allows all x-ray procedures even in challenging ER and outdoor environments. **MOVIX 4/8 DREAM** opens new horizons without compromising image quality.



3 STEPS READY MOBILE X-RAY UNIT



02



03



MOVIX 4/8 DREAM



TOUCHSCREEN CONTROL

Easy to use, the intuitive user interface allows programming of almost unlimited anatomical protocols. High-quality preview images are displayed **instantaneous** prior to transmission using the industry standard DICOM interface.

Advanced communication features have been added reducing the time required to share clinical information.

- USB devices
- CD/DVD Burner
- Ethernet cable
- Wifi

The integrated digital technology provides the highest image quality at the lowest dose.

A large range of Flat Panel Detectors is available to suit to your needs and provides optimum flexibility :

- **Multipurpose solution** with portable and detachable Flat Panel Detectors available in several sizes
- **Wireless solution** for a maximum of freedom with slim, lightweight Flat Panel Detectors.

ERGONOMICS

The dual side hand grips facilitate the tube positioning for efficient access to bedside and polytrauma patients.

Its large wheels* make **MOVIX 4/8 DREAM** fully dedicated to outdoor situations, operation fields, ER situations and sport events.

* option

With the **MOVIX 4/8 DREAM**, you enter in the digital era.

The wide fork* increases the stability of the **MOVIX 4/8 DREAM**. Its innovative design facilitates the accessibility where space is tight.



MOVIX 4/8 DREAM

IMAGE MANAGEMENT

Managing the explosion of digital data in healthcare is becoming increasingly complex. The need for an effective approach to information management and storage has never been greater and **STEPHANIX** has this.


STEPHANIX has far exceeded the basic functions of conventional storage to offer unique information management solutions with high-level features that focus on more efficient management and preservation of your data.

The **DReam View*** is an additional multimodality workstation, which provides diagnosis tools, CD/DVD burning functionality and DICOM 3.0 compliant capability.

* option



DX modalities

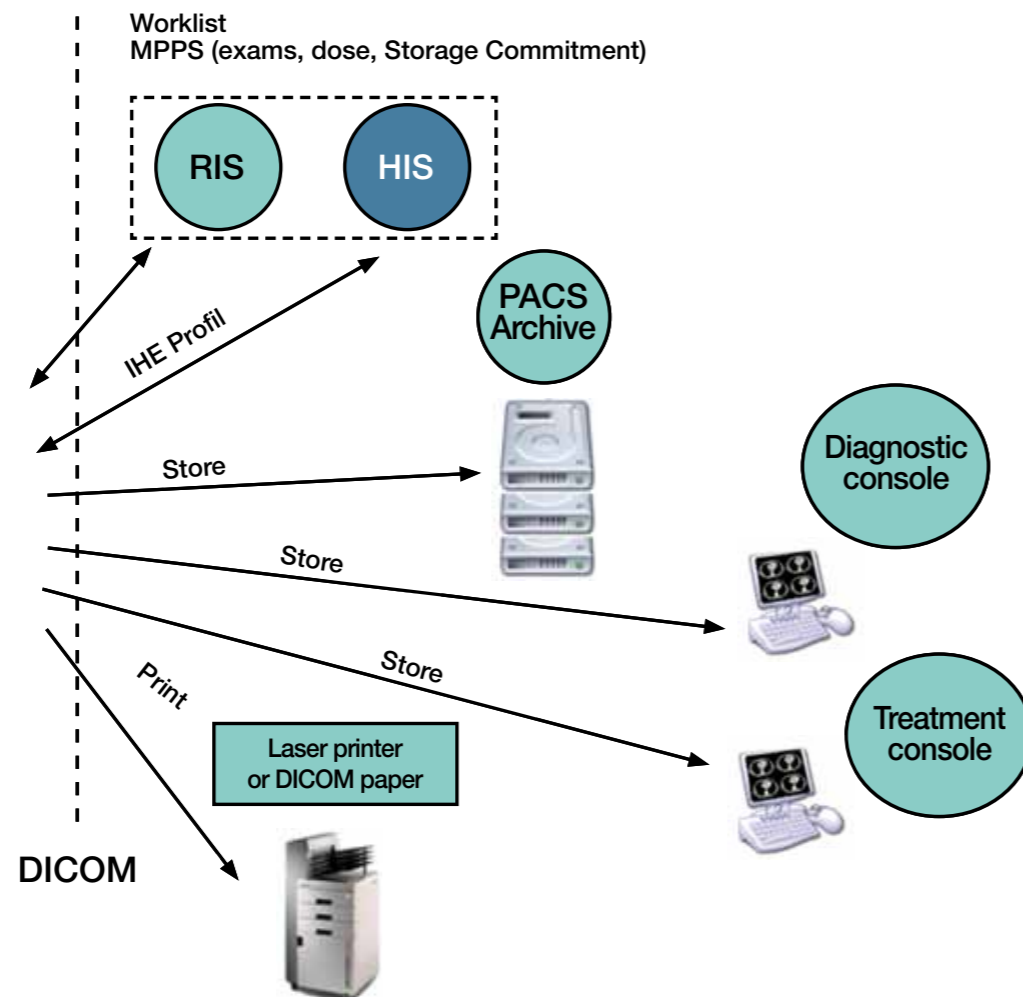

Example of the IT process.

APPLICATION

As safety comes first, we bring support to you to promote a special project with **Movix 4/8 DReam**. The design of the **Movix 4/8 DReam** is well adapted to respond to situations such as :

- **Military and Forensic**
Some cargo crates are available for delivery in air or by truck. Thanks to the robustness of the Flat Panel Detector, **MOVIX 4/8 DReam** can be used in hostile environment.
- **Outdoor and disaster**
Rapid and lightweight the **Movix 4/8 DReam** is ready for action with a simple 230 Vac to produce high quality imaging on the spot.



Latest generation digital radiography solutions with Integrated, Dynamic, Portable, Wired and Wireless detectors



STATIF PRO DREAM
Universal systems



Head Office
Rue Jean Moulin - ZI du Bayon
42150 La Ricamarie - FRANCE
Tel : + 33 4 77 47 81 60
contact@stephanix.com



www.stephanix.com

MOVIX 4/8 DREAM
Export department
15, rue du Saule Trapu
91882 Massy cedex - FRANCE
Tel : + 33 1 60 13 56 25
Fax : + 33 4 77 37 55 19