

MOVIX SERIES E+



STEPHANIX
MEDICAL IMAGING SOLUTIONS



Movix Series E + mobile units are capacitor assisted. They are available in three different powers: 16, 20 and 32 kW. Their different high kVp voltages enable exams in operating rooms, in emergency or intensive care units.



Thanks to its compactness, its large rear wheels and compact handle, **Movix Series E +** can move in small spaces and is even able to make a U-turn in less than 1m²




The design of the telescopic arm on a rotating column provides a better positioning, for exposures (end of bed for example), compared to a standard folding arm system.

The collimator movements on all axes with rotation of the column up to 300° together with the variable focal distance make positioning easy.



The soft keys console with 492 programmable APR controls the high frequency generator.

The programming consists in customised configurations with dedicated exposure parameters for each anatomical protocol.

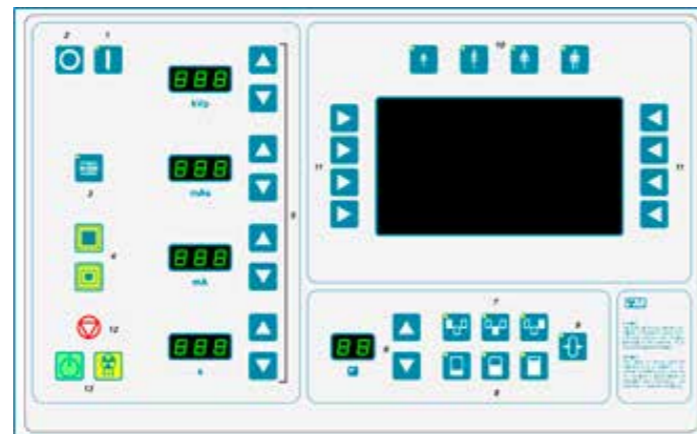

Handy and powerful (up to 32 kW), **Movix Series E +** is suitable for any kind of conventional procedures.

Max.
32kW

LATEST GENERATION OF DIGITAL RADIOGRAPHY SOLUTIONS: INTEGRATED, DYNAMIC, PORTABLE, WIRED AND WIRELESS DETECTORS.

D2RS
Dynamic Imaging



XTREME
General X-Ray Rooms



MOVIX SERIES DREAM
Mobile X-ray units



OMNISCOPE DREAM-S
Surgical C-arms



STEPHANIX
MEDICAL IMAGING SOLUTIONS

Head office : 10, rue Jean Moulin
Z.I. du Bayon - 42150 La Ricamarie - FRANCE
Tel. +33 4 77 47 81 60 - Fax: +33 4 77 37 55 19
export@stephanix.com

www.stephanix.com